# Final Project: Real-time electrocardiography (ECG)

#### Introduction:

The objective of this project is to measure heart rate in real time from subject electrocardiograph signal.

### ECG acquisition system:

ECG signals are collected from electrodes located on the right arm (RA), left arm (LA) and right leg (RL) of a subject (figure 1). The signals are amplified by an instrumentation amplifier and then filtered by a NOTCH filter (figure 2) before analog to digital conversion (ADC) using Teensy 3.2 development board. The ADC digitalized signals are then sent to the computer via serial port and displayed with a graphical user interface (GUI).

#### **Electrodes:**

Electrodes are placed on the right arm (RA), left arm (LA) and right leg (RL) of the subject. They are connected through wires to the input of the instrumentation amplifier, Vra and Vla, and the output of the driven right leg circuit, Vrl (figure 1).

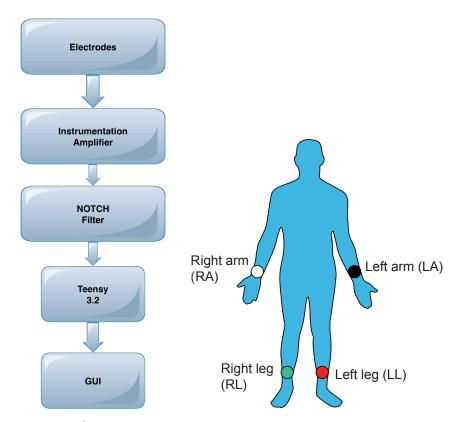


Figure 1. System diagram and electrode locations.

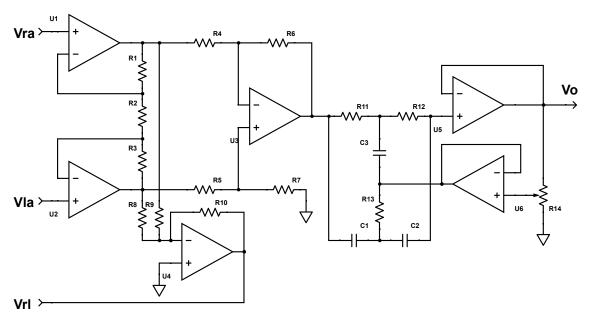


Figure 2. Instrumentation amplifier with driven right leg circuit and NOTCH filter.

# Instrumentation amplifier:

The amplifier gain has been set in 1156 (61dB). The resistor values used in the instrumentation amplifier and the right driven leg circuit are in table 1.

# Instrumentation amplifier

Value (Ohms)

39000

39000 1156

	Value (Ollilis)	
Resistor	±20%	
R1	39000	
R2	220	
R3	39000	
R4	12000	
R5	12000	

Right Driven Leg Circuit

	Value (kOhms)	
Resistor	±20%	
R8	39	
R9	39	
R10	220	

Table 1. Resistors used in instrumentation amplifier and right driven leg circuit.

#### **NOTCH filter:**

R6

R7

Gain

A second order NOTCH filter is used to diminish 60 Hz noise signal from power source. The element values are presented in table 2. The obtained cutoff frequency was 60.5 Hz.

NOTCH filter		
Cutoff Freq (Hz)	60.5	
R11 (Ohm)	5600	
R12 (Ohm)	5600	
R13 (Ohm)	2800	
C1 (uF)	0.47	
C2 (uF)	0.47	
C3 (uF)	0.94	

Table 2. NOTCH filter capacitance and resistor values.

#### ADC:

The ADC in the Teensy 3.2 board (figure 3) samples the analog signal at 1 kHz and sends the digital values to a computer via serial port. The Arduino code (appendix A) stores the ADC data in a buffer of size 'numofblocks' and sends each buffer to the computer via serial port.

### **Serial Port:**

The Teensy 3.2 serial port is set to work at 57,6 kbauds.



Figure 3. Teensy 3.2 development board.

#### GUI:

In the computer, a GUI has been implemented in Python 3.5 (figure 4). The GUI allows the visualization of the digitalized signal and its spectra (1<sup>st</sup> row), the filtered signal after a band-pass filter and it spectra (2<sup>nd</sup> row) and also the Pan and Tompkins peak finding algorithm and their locations (3<sup>rd</sup> row). The bottom row of the GUI has buttons and boxes that allow the user to start/stop the data reading, clear the buffered data, select low-pass and high-pass filter orders and cutoff frequencies, the window size in seconds to be visualized and the heart beat per minutes.

### Digital signal processing software:

The digital signals, sent by the hardware, are filtered to remove baseline and high frequency noise using programmable low-pass and high-pass filters (figure 4), in which the order and cutoff frequencies of the filters can be chosen by the user on the GUI. Then the signal is prepared to detect the R-peaks of the ECG by using the Pan and Tompkings' algorithm. The peak detection was performed using PeakUtils library.

The heart beat estimation is calculated by the inverse of the time distance between detected peaks multiplied by the ADC sampling rate and 60.



Figure 4. Digital signal processing implemented in Python.

## ECG reading using the system:

Figure 5 shows an example of an ECG signal measured by the implemented system. The top left window shows the digitalized ECG signal sent by the Teensy 3.2 board; on the left, the respective spectra is shown. The second row of windows display the signal after applying a band-pass filter on the raw data (left), in this case the low and high cutoff frequencies are 1 and 30 Hz respectively, with a low-pass filter order of 5 and high-pass filter order of 1. The respective spectra is displayed on the right column and it can be seen the low and high frequency attenuations. The third row displays on the left the Pan and Tompkins algorithm output after applying a derivative, signal squared and window integration, and on the right side the peak locations for the detected R-peaks of the ECG. The estimated hear beat is shown on the bottom-left part of the GUI, which for the example is 114 heart beats per minute.

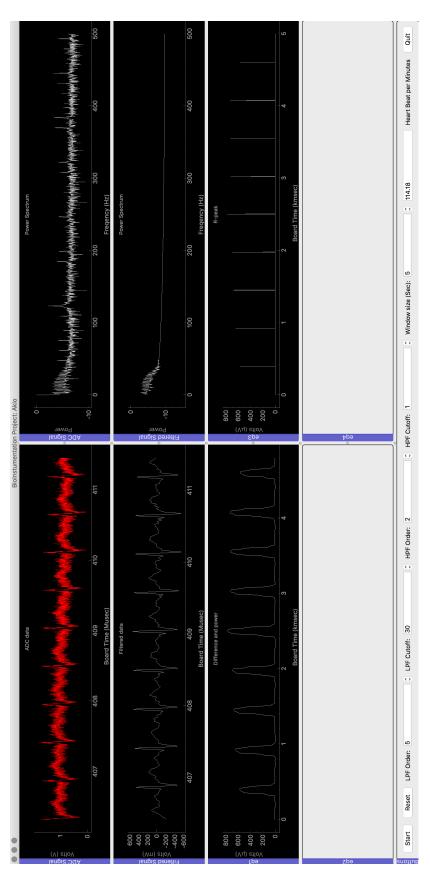


Figure 5. Graphical user interface in Python showing ECG signal, filtered signal and Pan and Tompkins algorithm output.

```
Appendix A: Arduino Code
// Arrays to save our results in
unsigned long values;
unsigned long time m;
int numofblocks = 200;
int delay const = 1000;
// Setup the serial port and pin 2
void setup() {
 // 300, 600, 1200, 2400, 4800, 9600, 14400,
 // 19200, 28800, 38400, 57600, or 115200
 Serial.begin(57600);
void loop() {
 for(int k=0;k<numofblocks;k++){</pre>
  time m = micros();
  values = analogRead(A0);
  Serial.print(time m);
  Serial.print(" ");
  Serial.print(values);
  Serial.print(" ");
  delayMicroseconds(delay const);
 Serial.print("\n");
Appendix B: Python Code
import numpy as np
import pygtgraph as pg
from pygtgraph.Qt import QtGui, QtCore
from pyqtgraph.dockarea import *
import serial
import scipy.signal as signal
import peakutils as pu
epsilon = np.finfo(float).eps
app = QtGui.QApplication([])
win = QtGui.QMainWindow()
# win = pg.GraphicsWindow(title="BioInstumentation Project: Akio")
area = DockArea()
win.setCentralWidget(area)
win.resize(1800, 950)
win.setWindowTitle('BioInstumentation Project: Akio')
# Create docks, place them into the window one at a time.
d11 = Dock("ADC Signal", size=(900, 300))
d12 = Dock("Filtered Signal", size=(900, 300))
```

```
d13 = Dock("Pan and Tompkins", size=(900, 300))
d15 = Dock("Buttons", size=(1800, 50))
d21 = Dock("ADC Signal", size=(900, 300))
d22 = Dock("Filtered Signal", size=(900, 300))
d23 = Dock("Peak Location", size=(900, 300))
area.addDock(d11, 'top')
area.addDock(d12, 'bottom', d11)
area.addDock(d13, 'bottom', d12)
area.addDock(d15, 'bottom', d13)
area.addDock(d21, 'right', d11)
area.addDock(d22, 'right', d12)
area.addDock(d23, 'right', d13)
# Widget plot p1:
p11 = pg.PlotWidget(title="ADC data")
p11.setLabel('left', "Volts", units='V')
p11.setLabel('bottom', "Board Time", units='usec')
p11.enableAutoRange('x', True)
p11.disableAutoRange('y')
p11.setYRange(0, 4)
curve11 = p11.plot(pen="r")
d11.addWidget(p11)
# Widget plot p3:
p21 = pg.PlotWidget(title="Power Spectrum")
p21.setLabel('left', "Power")
p21.setLabel('bottom', "Fregency", units='Hz')
p21.disableAutoRange('v')
p21.setYRange(-10, 0)
curve21 = p21.plot(pen="r")
d21.addWidget(p21)
# Widget plot p2:
p12 = pg.PlotWidget(title="Filtered data")
p12.setLabel('left', "Volts", units='V')
p12.setLabel('bottom', "Board Time", units='usec')
p12.enableAutoRange('x', True)
curve12 = p12.plot(pen="q")
d12.addWidget(p12)
# Widget plot p4:
p22 = pg.PlotWidget(title="Power Spectrum")
p22.setLabel('left', "Power")
p22.setLabel('bottom', "Frequency", units='Hz')
p22.disableAutoRange('v')
p22.setYRange(-15, 0)
curve22 = p22.plot(pen="q")
d22.addWidget(p22)
# Widget plot p2:
```

```
p13 = pg.PlotWidget(title="Difference, Squared and Integral")
p13.setLabel('left', "Volts", units='V')
p13.setLabel('bottom', "Board Time", units='msec')
p13.enableAutoRange('x', True)
curve13 = p13.plot(pen="c")
d13.addWidget(p13)
# Widget plot p2:
p23 = pg.PlotWidget(title="R-peak")
p23.setLabel('left', "Volts", units='V')
p23.setLabel('bottom', "Board Time", units='msec')
p23.enableAutoRange('x', True)
curve23 = p23.plot(pen="c")
d23.addWidget(p23)
# Widget for buttons:
p15 = pg.LayoutWidget()
B1 = QtGui.QPushButton('Start')
B2 = QtGui.QPushButton('Reset')
spinbox1 = pg.SpinBox(value=5, bounds=[1, 10], step=1)
spinbox12 = QtGui.QLabel("LPF Order:")
spinbox2 = pg.SpinBox(value=30, bounds=[5, 100], step=1)
spinbox22 = QtGui.QLabel("LPF Cutoff:")
spinbox3 = pq.SpinBox(value=2, bounds=[1, 6], step=1)
spinbox32 = QtGui.QLabel("HPF Order:")
spinbox4 = pg.SpinBox(value=1, bounds=[0.1, 10], step=0.1)
spinbox42 = QtGui.QLabel("HPF Cutoff:")
spinbox5 = pg.SpinBox(value=100, bounds=[10, 1000], step=1)
spinbox52 = QtGui.QLabel("Integral Window size:")
spinbox6 = pg.SpinBox(value=5, bounds=[2, 60], step=1)
spinbox62 = QtGui.QLabel("Window size (Sec):")
TBox1 = QtGui.QLineEdit()
TBox12 = QtGui.QLabel("Heart Beat per Minutes")
B100 = QtGui.QPushButton('Quit')
p15.addWidget(B1)
p15.addWidget(B2)
p15.addWidget(spinbox12)
p15.addWidget(spinbox1)
p15.addWidget(spinbox22)
p15.addWidget(spinbox2)
p15.addWidget(spinbox32)
p15.addWidget(spinbox3)
p15.addWidget(spinbox42)
p15.addWidget(spinbox4)
p15.addWidget(spinbox52)
p15.addWidget(spinbox5)
p15.addWidget(spinbox62)
p15.addWidget(spinbox6)
```

```
p15.addWidget(TBox1)
p15.addWidget(TBox12)
p15.addWidget(B100)
d15.addWidget(p15)
start stop = False
# Show Window
win.show()
# Serial port: OSX
raw = serial.Serial('/dev/cu.usbmodem1507501', 57600)
data0 = []
data1 = []
def update():
  global curve11, curve12, curve13, curve21, curve22, curve23
  global data0, data1
  numofblocks = 200 # Number of blocks in a package
  rate = 1000 # sampling rate in HZ
  window size = int(spinbox6.value()) # Window size in seconds
  integral win = int(spinbox5.value())
if start stop:
     xaxissize = window size * rate
     line = raw.readline()
     line val = [float(val) for val in line.split()]
     if len(line val) == 2 * numofblocks:
       for k in range(0, numofblocks):
          line val[2 * k] = line val[2 * k]
          line_val[2 * k + 1] = 3.3 * (line_val[2 * k + 1] - 1) / 1023
       if (not data0) | (len(data0) < xaxissize):
          for k in range(0, numofblocks):
             data0.append(int(line val[2 * k]))
             data1.append(line val[2 * k + 1])
       else:
          for k in range(0, numofblocks):
            data0.pop(0)
            data1.pop(0)
             data0.append(int(line val[2 * k]))
             data1.append(line val[2 * k + 1])
       xdata = np.array(data0, dtype='int')
       ydata = np.array(data1, dtype='float')
       curve11.setData(xdata, ydata)
       ydata filt = low pass filter(ydata, spinbox2.value(), rate,
spinbox1.value())
       ydata filt = high pass filter(ydata filt, spinbox4.value(), rate,
spinbox3.value())
       curve12.setData(xdata[integral_win:], ydata_filt[integral_win:])
       # ptr += 1
```

```
app.processEvents()
       if len(xdata) > 2*integral win:
          tdata = range(0, xaxissize)
          fx, fy = fouriertransform(tdata, ydata)
          curve21.setData(fx, fy)
          fx, fy = fouriertransform(tdata, ydata filt)
          curve22.setData(fx, fv)
          ydata_peak = find_peak(ydata_filt, integral_win)
          curve13.setData(ydata peak)
          ydata max = np.zeros(xaxissize)
          indexes = pu.indexes(ydata_peak, thres=0.5, min_dist=100)
          ydata max[indexes] = ydata peak[indexes]
          curve23.setData(tdata, ydata max)
          # print(np.diff(indexes))
          heart beat = 60*rate/np.mean(np.diff(indexes))
          TBox1.setText(str(np.round(heart_beat,decimals=2)))
timer = QtCore.QTimer()
timer.timeout.connect(update)
timer.start(0)
def find peak(x, n):
  y = np.diff(x)
  y = np.power(y, 2)
  y = np.cumsum(np.insert(y, 0, 0))
  y = (y[n:] - y[:-n]) / n
  y = y[100:]
  return y
def low pass filter(x, hf, rt, order):
  lowpass = signal.butter(order, hf/(rt/2.0), 'low')
  v = signal.filtfilt(*lowpass, x)
  return y
def high pass filter(x, lf, rt, order):
  highpass = signal.butter(order, lf/(rt/2.0), 'high')
  y = signal.filtfilt(*highpass, x)
  return y
def fouriertransform(x, y):
  # Perform fourier transform. If x values are not sampled uniformly,
  # then use np.interp to resample before taking fft.
  dx = np.diff(x)
  uniform = not np.any(np.abs(dx - dx[0]) > (abs(dx[0]) / 1000.))
  if not uniform:
     x2 = np.linspace(x[0], x[-1], len(x))
```

```
y = np.interp(x2, x, y)
     x = x2
  f = np.abs(np.fft.fft(y) / len(y)) ** 2
  y = np.log10(f[0:len(f) / 2] + epsilon)
  dt = (x[-1] - x[0]) / 1000
  x = np.linspace(0, 0.5 * len(x) / dt, len(y))
  return x, y
def start_scope():
  global start stop
  if start stop:
     start stop = False
     B1.setText("Start")
  else:
     start stop = True
     B1.setText("Stop")
B1.clicked.connect(start scope)
def reset scope():
  global data0, data1
  data0 = []
  data1 = []
B2.clicked.connect(reset scope)
spinbox5.valueChanged.connect(reset scope)
def close scope():
  raw.close()
  exit()
B100.clicked.connect(close scope)
def clear data():
  global xdata, ydata
  xdata = []
  ydata = []
spinbox5.sigValueChanging.connect(clear_data)
if __name__ == '_ main ':
  import sys
  if (sys.flags.interactive != 1) or not hasattr(QtCore, 'PYQT VERSION'):
     QtGui.QApplication.instance().exec ()
```